

CARENA

For Day Surgery

Mobile like a stretcher,
comfortable as a bed

DAY SURGERY



Mobile like a stretcher, comfortable as a bed

Specially designed for day hospital needs, the 4-sectional Carena Day Surgery stretcher provides comfort for both patients and clinical staff.

Great mobility and patient safety


Easy to move around also in narrow corridors, the stretcher is ideal for transporting patients/clients safely and conveniently from the operating room to recovery.

Easy to clean and maintain in good condition

The mattress base plates are easy to wash and remove. Thanks to the structure of the bed, maintenance of components is made easy, too.

Equipped for day hospital needs

Naturally, the stretcher is equipped with all necessary accessories that you'll need on your day hospital stretcher, such as convenient side rails, infusion rod, monitor holder, journal basket and lifting pole.

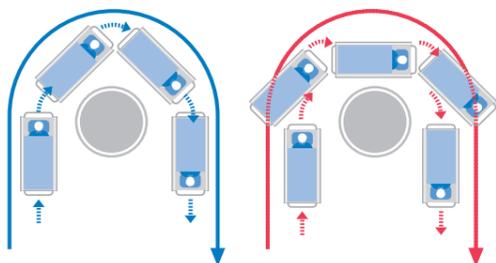

Improved mobility

The turning radius decreases with the help of the fifth wheel (picture left) thus increasing mobility. Without the wheel the radius is larger (picture right).



No need for separate recovery furniture

Long and wide enough, and with adjustable seating, Carena Day Surgery doubles as a recovery bed/chair providing great comfort also for the patient.

During recovery, day surgery clients can adjust the seat themselves with the hand control to find the most comfortable position before leaving the clinic.

As it goes really low, also shorter clients can easily enter and leave the bed once they are ready to take a walk. They'll love it!



Technical features:

- Takes patients up to 200 kg
- Safe working load (SWL) is 230 kg
- Minimum height 390 mm, maximum height 900 mm
- Width 800 mm
- Electricity-powered height and mattress base adjustment
- Both electricity-powered and manual Trendelenburg adjustment: -12° to +6°
- Zero level position in hand control
- Fifth wheel for improved mobility
- Durable construction, chrome base and white bed ends
- The Carena Day Surgery stretcher meets all the newest directives and regulations for medical devices: EN 60601-1, EN 60601-1-2, EN 60601-2-52.



"We really like the Carena. It is so easy to manoeuvre and it fits well into our new day hospital. Additionally, also tall patients can feel comfortable and safe on it, which is a great plus!"

Katleen Devriendt, Head Nurse at Sint Blasius Day Hospital

Easy to manoeuvre

Safe for patients but small enough for transportation

Comfort for clients in recovery


For further details, contact your nearest Merivaara representative, see www.merivaara.com for contact info.

Solutions for efficient physical patient flow and well-being

Merivaara provides a wide range of hospital-grade furniture such as operating tables, medical lights, examination tables, trolleys or stretchers for transportation and day surgery, as well as birthing and patient beds. Our versatile product offering enables convenient and cost-effective physical patient flow within the hospital, in day surgery and health clinics, as well as in nursing homes and home care. Our priority is to enhance comfort and well-being of both patients and their caregivers.

Established in 1901, Merivaara has more than 100 years of experience in designing and manufacturing hospital-grade furniture. Today, Merivaara's patient beds, surgical tables and trolleys are highly appreciated by users in more than 100 countries thanks to their ease-of-use, durability and ergonomic design.

Merivaara's products are CE-labelled and the company's quality management system complies with EU directives for medical device, and is certified by DNV according to ISO 9001:2008 and ISO 13485:2003 standards. The company's environmental management system is ISO 14001:2004 certified.

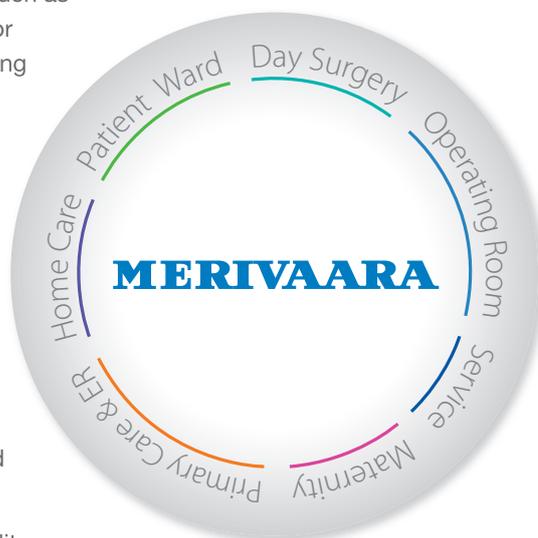

Merivaara Corp.

Puustellintie 2, FI-15150 LAHTI, FINLAND
Tel. +358 3 3394 611, fax +358 3 3394 6144
merivaara@merivaara.fi
www.merivaara.com

The status of the information, specifications and illustrations in this publication is indicated by the date of publication. Whilst the greatest care has been taken over its compilation, Merivaara Corp. is not responsible for any errors or omissions. Merivaara Corp. reserves the right to make changes to the technology, features, specifications and design of the equipment and models without notice. Please check the Merivaara website for the latest updated specifications and other information.

All trademarks are Merivaara property unless otherwise stated. Patents granted or pending.