

Therapy Chair Série II



ActualWay® Therapy Chair - Série II
Reliability, Stability and Performance
This Medical Treatment Chair has been developed to give comfort, security and wellness during long treatments such as Dialysis, Chemotherapy, Oncology and Blood Transfusions

FEATURES AND ADVANTAGES

- Easy to use and maintain
- Excellent delivery time and after sales service
- 3 or 4 motors
- Adjustable head cushion
- Handset including the Trendelenburg operation key
- Arm rest vertical and horizontal position adjustment
- Arm rest inclination adjustment
- Built-in emergency/crash button and very easy to activate to get the Trendelenburg position -12° in 12 seconds
- Granted Stability
- Tested and approved load capacity: 200Kg + the optional accessories
- Central locking system on 10cm diameter castors
- Rotating castors with medical properties (Polypropylene, thermoplastic rubber, grey non-marking, precision ball bearing, solid plug fitting)
- High quality electrical components and accessories
- Epoxy painting metal frame
- Included battery
- Special fireproof, waterproof and efficient against bacteria (MRSA) French coated fabric and easy to clean, class M2
- Range of 15 available colours

LARGE RANGE OF OPTIONS

- Weight Scale M and Class III
- Optional Special Dimensions on the 4 motors version Chair: Length 217cm and width 60cm
- Extra comfort upholstery
- Special density foam
- IV serum pole
- Foldable and adjustable, 5Kg capacity, rectangular tray in ABS
- Paper roll holder
- Back handle
- Seamless upholstery
- Upholstery colour range: class M1 (Urban, Sirius,...)
- ABS protection bottom cover
- 1 110V control box
- Foot rest motor
- DVD player
- Reading Lamp

STANDARD DIMENSIONS


- 205x54cm
- Total width (including the arm rest): 88cm
- Frame: 105x65x25cm
- Height: 55-75cm

TECHNICAL SPECIFICATIONS

- Power Supply: 230V /50-60Hz
- Protection Category: IP44
- Handset protection category: IP66

LOGISTIC DATA

- Net weight: 85Kg
- Palett: 0,80x1,20x1,60m
- Gross weight: 120Kg

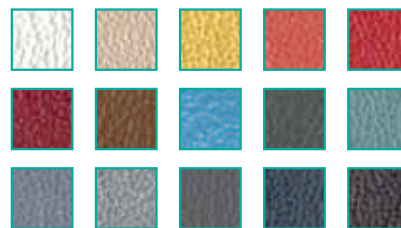

Distributor:



Producer:
VIA ACTUAL-COMÉRCIO INTERNACIONAL, LDA
CE Castelo da Maia
Rua Manuel Assunção Falcão, 481- Fracção H
4475-041 Santa Maria de Avioso - Maia
PORTUGAL

Tlf. / Phone: +351 229 829 663/4
Fax. +351 229 829 665
E-mail: comercial@actualway.pt
comercial@actualway.eu
www.actualway.pt